


SafeGuardPro™

Ductless Fume Cupboard

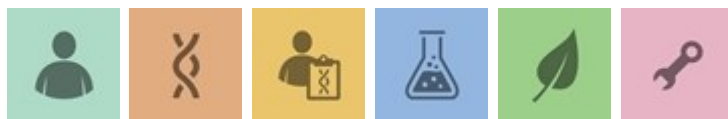

The SafeGuardPro™ ductless fume cupboard provides both operator and environment protection for corrosive and solvent-based applications. With a portable, modern design made here in Australia, this polypropylene fume cabinet offers exceptional chemical resistance—a perfect alternative to standard metal fume cupboards which risk rust and reduced safety. Minimal installation is required, and with no ductwork, the SafeGuardPro™ is an economical, efficient option for your lab.

A state-of-the-art control system provides airflow and filter monitoring, ensuring user protection 24 hours a day. Multiple filtration options, including HEPA and carbon filters, are available—our specialists can help you in determining the best solution for your application. The SafeGuardPro™ ductless fume cabinet exceeds requirements of AS/NZS 2243.9:2009 Safety in laboratories – Recirculating fume cabinets.





Features

- Ductless, portable design
- Designed & made in Australia
- Wide range of filtration options available
- Minimal installation – will plug directly into standard 220V 10AMP electrical outlet
- Toughened glass sliding sash
- Polypropylene construction
- Features automatic safety controller
- Low airflow & filter replacement alarm

SafeGuardPro™

Ductless Fume Cupboard

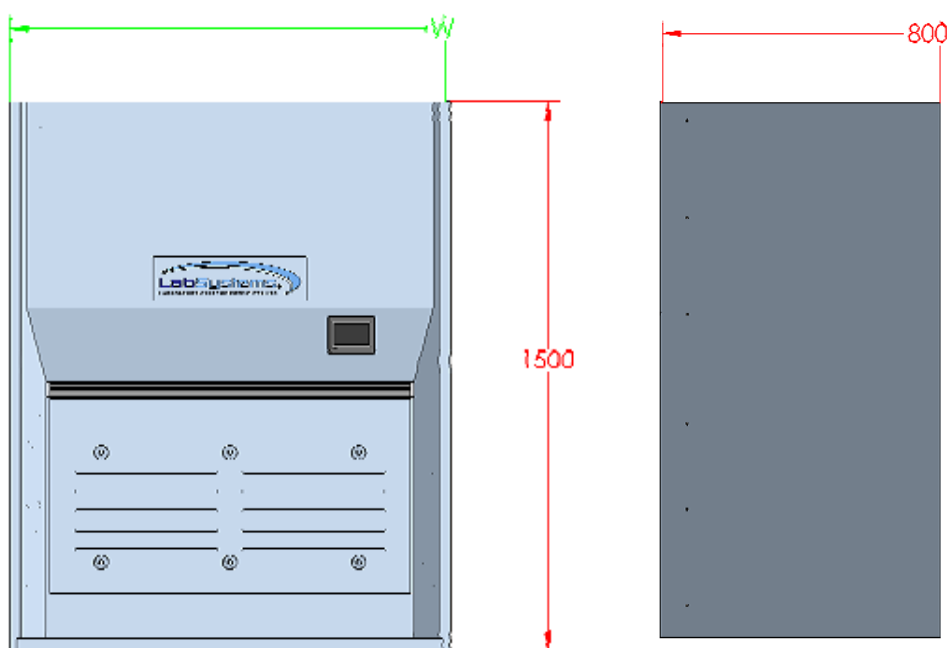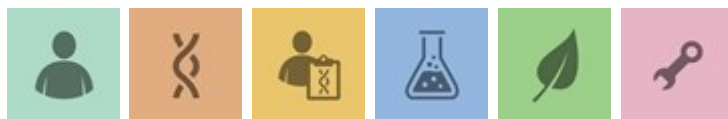

Optional Features

- Polypropylene sink & drain
- Gas & water ports
- Cord port
- Custom base inserts
- Steel under bench frame
- Double filter application
- Electric height-adjustable bench

DIMENSIONS	SG Pro900	SG Pro 1200	SG Pro 1500	SG Pro 1800	SG Pro 2000	SG Pro 2400
HEIGHT	1500mm	1500mm	1500mm	1500mm	1500mm	1500mm
WIDTH	900mm	1200mm	1500mm	1800mm	2000mm	2400mm
DEPTH	800mm	800mm	800mm	800mm	800mm	800mm